AT51140/AT51141
INSTRUKCJA UŻYTKOWNIA, LASKA
INSTRUCTIONS FOR USE, CANE
GEBRAUCHSANWEISUNG, GEHSTOCK
NÁVOD K POUŽITÍ, HŮL
NÁVOD NA POUŽITIE, PALICA
NOTICE D'UTILISATION, CANNE
GEBRUIKSAANWIJZING, WANDELSTOK
INSTRUCCIONES DE USO, BASTÓN
ISTRUZIONI PER L'USO, BASTONE
BRUKSANVISNING, KÄPP

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En

Thank you for purchasing our product. Please read these instructions for use carefully. The product must not be used before reading and understanding this manual. If you do not understand the warnings, notes, or recommendations, contact a healthcare professional or the seller to avoid damage to the product or personal injury. **WARNING:** All parts should be checked for damage that may have occurred during transport. If such damage is noticed, the product **MUST NOT** be used. For more information, contact the seller.

TARGET PATIENT GROUP

Persons suffering from diseases, dysfunctions, or injuries for which this device is intended for treatment, rehabilitation, or compensation (see the section concerning the intended use of the device in this manual). The product may be purchased independently by the user or based on the recommendation of a physician, therapist, or another specialist. Whether the product is purchased independently or recommended by a healthcare professional, the available sizes, required functions, product variants, indications and contraindications for use, and the information provided by the manufacturer should be taken into account.

INTENDED USE: The product is intended to provide support while walking for persons with reduced mobility.

INDICATIONS

The product may be used in particular in the following cases:

balance and stability disorders while walking,
reduced muscle strength in the lower limbs,
recovery period after injuries or surgical procedures of the lower limbs,
musculoskeletal disorders that make independent movement difficult,
neurological conditions affecting gait and balance,
elderly persons requiring additional support while walking.

CONTRAINDICATIONS

physical or mental limitations (e.g., visual impairment) that prevent safe use of the product,
inability to maintain balance independently while standing,
severe coordination disorders preventing safe use of the cane,
upper limb disorders preventing safe support on the cane,
the need to use more advanced mobility aids (e.g., rollator, walker) due to severe mobility limitations,
skin damage, open wounds, or severe pain at the hand support area that prevents the use of the cane.

LIST OF COMPONENTS: Cane, instruction manual

PRODUCT FEATURES

AT51140 – FOLDABLE CANE (four segments)

AT51141 – STANDARD CANE

Lightweight and durable aluminum construction

Height adjustment allowing the cane to be adapted to the user's height

Ergonomic handle ensuring comfortable and stable grip

Anti-slip rubber tip improving safety and stability during walking

Designed for everyday mobility support

Corrosion-resistant materials suitable for long-term use

TECHNICAL SPECIFICATION AT51140

Height adjustment: 76–86 cm

Tube diameter (upper section): Ø22 mm

Tube diameter (lower section): Ø19 mm

Weight: 0.37 kg

Frame material: anodized aluminum

Handle material: PA (polyamide)

Tip material: anti-slip rubber

Maximum user weight: 100 kg

TECHNICAL SPECIFICATION AT51141

Height adjustment: 69–92 cm

Tube diameter (upper section): Ø22 mm

Tube diameter (lower section): Ø19 mm

Weight: 0.35 kg

Frame material: anodized aluminum

Handle material: PA (polyamide)

Tip material: anti-slip rubber

Maximum user weight: 100 kg

NOTE: Individual product dimensions may vary within ± 2 cm.

This symbol indicates the maximum user weight.

INSTRUCTIONS FOR USE

Adjusting the cane height

Wear the footwear in which the cane will most often be used and stand in a natural position.

Place the cane approximately 20 cm in front of the foot.


Adjust the cane height so that the handle is at hip joint level.

When holding the cane, the elbow should be slightly bent – approximately 30°.

To adjust the height, press the locking button, set the appropriate cane length, and ensure that the button locks properly in the adjustment hole.

Walking with the cane

The cane should be held on the side opposite the weaker or injured leg. While walking, move the cane forward together with the weaker leg. Then take a step with the stronger leg. Try to walk naturally, taking even and steady steps.

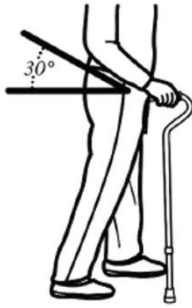**WARNING**

In the event of a serious incident related to the device that directly or indirectly led, could have led, or may lead to any of the following:

- a) death of a patient, user, or other person
- b) temporary or permanent serious deterioration of the health of a patient, user, or other person
- c) serious public health threat

the incident must be reported to the manufacturer and the competent authority of the Member State in which the user or patient resides.

WARNING: If pain, allergic reactions, or other concerning or unclear symptoms related to the use of the medical device occur, consult a healthcare professional.

WARNING: The product must not be used for purposes other than its intended use.

CLEANING AND MAINTENANCE

Clean the cane handle using mild soap, a gentle detergent, or a household cleaning agent. After cleaning, wipe the surface with a damp cloth and dry it.

Do not use solvents or aggressive chemicals such as acetone, paint thinner, lacquer thinner, or turpentine, as they may damage the product surface.

Before each use and after use, check the technical condition of the cane, especially:

- whether the cane structure is deformed or cracked,
- whether the handle is stable and not loose,
- whether anti-slip elements are in the correct position,
- whether the rubber tip is worn or damaged,
- whether the height adjustment button works properly and is not damaged.

If any damage is detected, the product must not be used until the defect has been repaired or damaged parts replaced.

STORAGE: Store the product in a dry place. Do not expose it to direct sunlight.

DISPOSAL: After the end of its service life, the product should be disposed of as normal household waste. Follow local waste segregation regulations. Protect the environment.

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Model:
Serial number:

Point of sale stamp

and a legible signature of the Seller:

1. ANTAR MEDICAL Sp. z o.o. ,03-068 Warsaw, ul. Zawisłańska 43, hereby grants a 12-month warranty on the product, starting from the date of purchase by the Buyer.
2. During the warranty period ANTAR undertakes to carry out the repairs necessary to restore proper functioning of the product, within 14 days from the date of receipt of the product by the service centre, together with the proof of purchase and the warranty card.
3. ANTAR reserves the right to replace the product if the repair is deemed uneconomic. Repair or replacement of faulty products does not extend the warranty period.
4. If the product is found to be malfunctioning, please contact your retailer/distributor or (if it is impossible) send it to ANTAR (as specified in point 1).
5. Before use, read the instructions carefully and follow them. If the product is used improperly and contrary to recommendations, the warranty will be irrevocably lost.
6. The warranty does not cover damage caused by external mechanical force, load exceeding the maximum allowed weight of the user, as described in the manual, as well as contact of the product with alcohol, greasy substances or petrol. The warranty does not cover parts which are subject to normal wear or tear during use, such as the bearing, rubber parts and others.
7. The product delivered to the service in connection with the complaint cannot cause an epidemic hazard.
8. THIS WARRANTY IS ONLY VALID TOGETHER WITH THE PROOF OF PURCHASE (BILL, RECEIPT, INVOICE). IN CASE OF A COMPLAINT, A PROOF OF PURCHASE MUST BE ATTACHED.
9. Complaints submitted without the proof of purchase and a warranty card with the product serial number will not be considered.
10. The warranty does not cover parts and components subject to natural wear and tear, such as wheels, tires, rubber pads, fuses, etc., the wear of which results from normal use of the product. The warranty period for batteries is 6 months.

ANTAR MEDICAL Sp. z o.o.
03-068 Warsaw, ul. Zawisłańska 43
Tel. 22 518 36 00, Fax 22 518 36 30
www.antar.net

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