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EN

Thank you for purchasing our product. Please read this manual carefully before use. Do not use the product until you have read and understood this manual. If you do not understand the warnings, notes, or recommendations, contact a healthcare professional or the seller to avoid damage to the product or injury.

WARNING:

Check all parts of the product for damage caused during transport. If any damage is found, DO NOT use the product. For more information, contact the manufacturer. Before putting on the product, ensure it is in the correct position, with the slots for the stays at the top.

APPLICATION AT04502/3/10 PRO

The spine brace improves lumbar spine stabilization, reduces pain, and enhances daily comfort. Stabilization is provided by flexible fastening straps and contoured reinforcements (stays), which relieve the tissues around the spine by redistributing the load forces on the spine structures. This prevents harmful changes in the lower spine. Thanks to its anatomical design, the brace fits the body perfectly and can be worn discreetly under clothing. High-quality Velcro fasteners ensure easy and precise adjustment.

APPLICATION AT04505 PRO

The spine brace improves stabilization of both the lumbar and thoracic spine, reducing pain and enhancing daily comfort. Stabilization is ensured by flexible fastening straps, contoured reinforcements, and straight stays, which relieve the tissues around the spine by redistributing load forces on spinal structures. Its anatomical design allows a perfect fit to the body, and it can be worn discreetly under clothing. High-quality Velcro fasteners ensure easy and precise adjustment.

INDICATIONS AT04502/3/5/10 PRO

Spinal disorders
Multilevel discopathy
Osteoporosis
Degenerative changes
Spinal instability
Lumbar spondylolisthesis
Compression fractures due to osteoporosis
Sciatica
Pain syndromes in the lumbar-sacral region
Post-operative rehabilitation

CONTRAINDICATIONS

Physical or mental limitations (e.g., vision impairment) preventing safe handling of the product. Do not use the product after applying warming gels or ointments.

FEATURES AT04502 PRO

Equipped with 6 contoured steel reinforcements: 4 (28 cm long) and 2 (25 cm long)
Features a handle system for easy application and removable stays for easy cleaning

FEATURES AT04503 PRO

Equipped with 4 contoured steel reinforcements, each 19 cm long
Features a handle system for easy application and removable stays for easy cleaning

FEATURES AT04510 PRO

Equipped with 4 contoured steel reinforcements, each 28 cm long
Features a handle system for easy application and removable stays for easy cleaning

FEATURES AT04505 PRO

Equipped with 4 contoured steel reinforcements (28 cm long) and 2 straight reinforcements (45 cm long)
Features a handle system for easy application and removable stays for easy cleaning

SIZES AT04502/3/5/10 PRO (*Measure waist circumference to determine the appropriate size*):

S: 68–78 cm
M: 76–86 cm
L: 82–92 cm
XL: 90–100 cm
XXL: 98–110 cm
XXXL: 108–125 cm

HEIGHT:

AT04502 PRO – 32 cm
AT04503 PRO – 25 cm
AT04510 PRO – 32 cm
AT04505 PRO – 55 cm

PACKAGE CONTENTS

Spine brace, user manual

WARNING:

In the event of a “serious incident” related to the product that has directly or indirectly caused, could have caused, or is likely to cause:

- a) Death of a patient, user, or another person;
- b) Temporary or permanent deterioration of health;
- c) Serious public health threat;

this “serious incident” must be reported to the manufacturer and the relevant authority of the member state where the user or patient resides. In Poland, the competent authority is the Office for Registration of Medicinal Products, Medical Devices, and Biocidal Products.

WARNING:

If pain, allergic reactions, or other concerning symptoms related to product use occur, consult a healthcare professional immediately.

WARNING:

It is prohibited to use the product for purposes other than those intended.

WARNING:

The manufacturer is not responsible for damage caused by neglecting cleaning and maintenance guidelines or failure to follow the recommendations in this manual.

INTENDED USER GROUP


This product is intended for individuals suffering from conditions, dysfunctions, or injuries requiring treatment, rehabilitation, or compensation (see the purpose section in this manual). The product can be purchased independently or based on a recommendation from a doctor, therapist, or specialist. Regardless of the purchase method, available sizes, required functions, product variations, indications, and contraindications should be considered, as well as information provided by the manufacturer.

USAGE AT04502/3/10 PRO

Ensure the brace is in the correct position. Place your thumbs in the belt pockets. Fasten the main strap tightly with Velcro, compressing the abdominal belt appropriately.

Fasten the two additional tensioning straps by pulling both ends simultaneously.

Finally, make sure the product fits properly to the spine, is positioned at the correct height, and the straps are tightened to provide adequate stabilization.

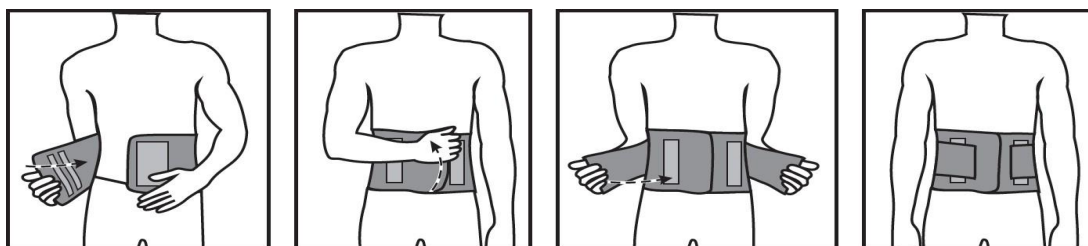
**USAGE AT04505 PRO**

Ensure the brace is in the correct position, with long straps over the shoulders. Place your thumbs in the belt pockets. Fasten the main strap tightly with Velcro, compressing the abdominal belt.

Fasten the two additional tensioning straps by pulling both ends simultaneously. (AT04510 PRO)

Pass the straps under the arms, cross them behind the back, pull them to the front, and fasten them on the main belt.

Finally, verify the brace properly fits the spine, is at the correct height, and the straps are securely tightened to provide stabilization.

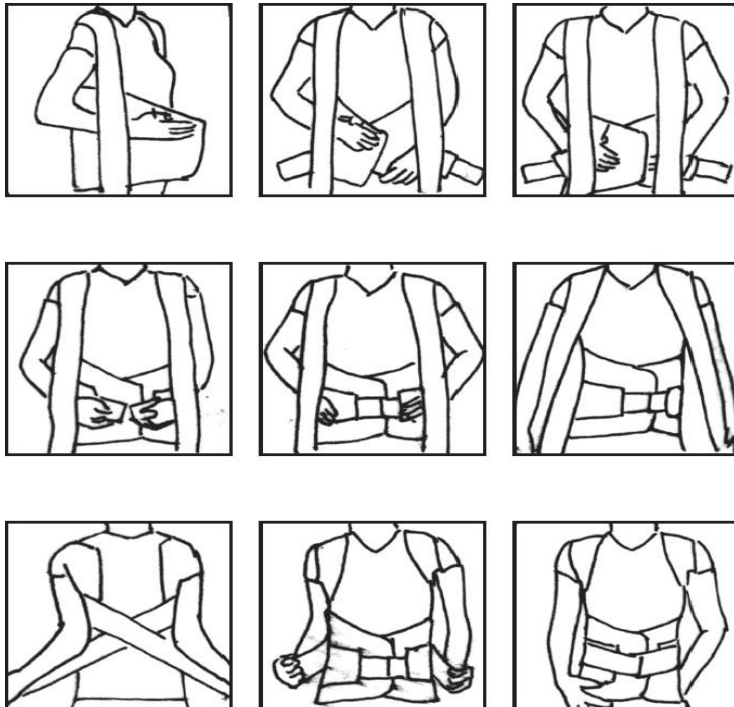CLEANING

Hand wash the product with mild detergent at a temperature of up to 40°C. Do not machine wash. Do not use strong detergents, bleach, or abrasive materials. Dry flat in the shade, avoiding direct sunlight. Do not tumble dry, dry clean, or iron.

STORAGE

Store the product in a dry, cool place, away from moisture and direct sunlight.

DISPOSAL

The medical device can be disposed of as regular household waste.

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Model:
Serial number:

Point of sale stamp

and a legible signature of the Seller:

1. ANTAR Sp. J., 03-068 Warsaw, ul. Zawiślańska 43, hereby grants a 12-month warranty on the product, starting from the date of purchase by the Buyer.
2. During the warranty period ANTAR undertakes to carry out the repairs necessary to restore proper functioning of the product, within 14 days from the date of receipt of the product by the service centre, together with the proof of purchase and the warranty card.
3. ANTAR reserves the right to replace the product if the repair is deemed uneconomic. Repair or replacement of faulty products does not extend the warranty period.
4. If the product is found to be malfunctioning, please contact your retailer/distributor or (if it is impossible) send it to ANTAR (as specified in point 1).
5. Before use, read the instructions carefully and follow them. If the product is used improperly and contrary to recommendations, the warranty will be irrevocably lost.
6. The warranty does not cover damage caused by external mechanical force, load exceeding the maximum allowed weight of the user, as described in the manual, as well as contact of the product with alcohol, greasy substances or petrol. The warranty does not cover parts which are subject to normal wear or tear during use, such as the bearing, rubber parts and others.
7. The product delivered to the service in connection with the complaint cannot cause an epidemic hazard.
8. THIS WARRANTY IS ONLY VALID TOGETHER WITH THE PROOF OF PURCHASE (BILL, RECEIPT, INVOICE). IN CASE OF A COMPLAINT, A PROOF OF PURCHASE MUST BE ATTACHED.
9. Complaints submitted without the proof of purchase and a warranty card with the product serial number will not be considered.

ANTAR Sp. J.
03-068 Warsaw, ul. Zawiślańska 43
Tel. 22 518 36 00, Fax 22 518 36 30
www.antar.net

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