AT51211

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En

Thank you for purchasing our product. Please read the instructions for use carefully. Do not use the device before reading and understanding these instructions. If you do not understand the warnings, notes and recommendations, contact a healthcare professional or the seller to avoid damage to the product or bodily injury. NOTE: Check all parts for damage. Before first use, thoroughly inspect the technical condition of the product and make sure that all parts are complete and correctly installed (if applicable, pay special attention to screws, nuts, rivets and other fastening and structural components). Make sure that the product is not defective or damaged. If any damage or missing parts are noticed, DO NOT use the product. For additional information, contact the seller.

TARGET PATIENT GROUP

Persons dealing with diseases, dysfunctions or injuries, for the treatment/rehabilitation or compensation of which this device is intended (see the section concerning the intended use of the device in these instructions for use). The product may be purchased independently by the user or based on the recommendations of a doctor, therapist or another specialist. Both in the case of independent purchase of the product and when indicated by a doctor/therapist/other specialist, the available sizes/necessary functions and variants of the device, indications and contraindications for use, as well as the information provided by the manufacturer, must be taken into account.

INTENDED USE: The electric raised toilet seat is intended to assist persons with reduced mobility when sitting down on and standing up from the toilet through controlled lifting and lowering of the user.

INDICATIONS:

Indications for use include, among others:

- limited mobility resulting from age,
- weakness of the muscles of the lower limbs,
- disorders of balance and body stability,
- orthopaedic problems making sitting down and standing up difficult,
- recovery period after injuries, surgical procedures or hospitalisation,
- neurological diseases affecting mobility and movement control,
- degenerative diseases of the hip joints, knee joints or spine,
- increased risk of falls when using the toilet,
- the need to reduce physical effort when performing everyday hygiene activities,
- the need for additional support and safety when using the toilet.

The product may be used in home conditions and care facilities in accordance with its intended use.

CONTRAINDICATIONS:

Physical or mental limitations (e.g. visual impairment) that prevent safe handling of the product.

Do not use the product:

- in persons with a body weight exceeding 150 kg,
- if the person is unable to independently maintain a sitting position without appropriate support,
- in persons with serious disturbances of consciousness or orientation preventing safe use of the device,
- in the case of unstable injuries or conditions requiring individual medical equipment,
- in the case of significant body deformities preventing proper use of the product,
- in the case of product damage or malfunction of the lifting mechanism,
- in a wet environment or one exposed to direct contact with water,
- by children, as a toy or climbing device,
- in a manner inconsistent with the manufacturer's intended use.

In the case of serious medical conditions or doubts regarding the possibility of using the product, consultation with a doctor or medical personnel is recommended.

GENERAL SAFETY WARNINGS

- Before use, assess whether the user is suitable to use the product, taking into account factors such as body weight, body build, cognitive abilities and mobility limitations. Users who are not suitable should not use the AX12 in order to avoid injury. Persons with disabilities should use the product only under the supervision of a caregiver or with the use of an alarm system.

- Before use, make sure that the ambient conditions meet the following requirements:

Operating temperature and humidity: 0°C to +45°C; relative humidity from 10% RH to 75% RH

Storage temperature and humidity: -10°C to +60°C; relative humidity from 0% RH to 75% RH

- Operating altitude: the controller is intended for operation at an altitude of up to 2000 metres. (At high altitudes)
- The packaging contains a plastic bag. To avoid the risk of suffocation, keep it out of the reach of children. The plastic bag is not a toy and should not be given to infants or small children.

WARNINGS CONCERNING USE

- The maximum load of the product is 150 kg. Exceeding this limit may damage the product or pose a risk to the user's safety.
- The product is intended for sitting only. Do not use it as a platform, ladder, support or for any other purposes inconsistent with the intended use of the product.
- During use, do not insert any objects, such as sticks, into moving parts. This may damage the product or create a safety hazard.
- When sitting down, do not apply excessive force to the seat. The user should sit down calmly and slowly, aiming for the centre of the seat. One hand should hold the handle and the other should operate the controller for lifting or lowering the product.
- Do not stand or kneel on the product during use.
- Do not disassemble or modify the product yourself, as this may damage it and create a safety hazard.
- After assembly, before proper use, perform an operating test to ensure that the raised toilet seat lifts and lowers correctly and that it has sufficient working space.
- Do not use small objects, such as brushes, which could be pulled into the mechanism. Do not place objects, such as clothing, near the seat of the raised toilet seat.
- Use with wet hands is strictly prohibited due to a potential safety hazard.

- Keep animals away from the product and from the device while it is in use to avoid injuries or other safety hazards.
- The raised toilet seat is not a toy. Do not allow children to operate or use the device.

WARNINGS CONCERNING BATTERY SAFETY

- The hand controller of the raised toilet seat must not be disassembled. Disassembly may damage the product or cause bodily injury.
- Only the hand controller supplied with the product should be used to lift and lower the device. Use of another controller may cause malfunction of the device or create a safety hazard.
- Do not place the product near heat sources to avoid potential damage or safety hazards.
- During use, make sure that the product's power plug remains dry and clean.
- Never operate the device with wet hands in order to prevent the risk of electric shock.

DECLARATION CONCERNING ELECTROMAGNETIC INTERFERENCE AND ELECTROMAGNETIC COMPATIBILITY (EMI/EMC)

The product has been designed and tested in accordance with the relevant international and national standards in order to ensure compliance with the requirements concerning electromagnetic interference (EMI) and electromagnetic compatibility (EMC). The details of the declaration are presented below:

1. General Declaration of Conformity

The product meets the requirements of the following standards concerning electromagnetic interference and electromagnetic compatibility:

EN 55014-1, EN 55014-2 (standards concerning the electromagnetic compatibility of household appliances, electric tools and similar devices).

2. User Information and Warnings

Using incompatible accessories or making unauthorised modifications may reduce electromagnetic compatibility and affect the operation of the device.

If the device is exposed to external electromagnetic interference or interferes with the proper operation of other devices, immediately change the position of the device or stop using it.

3. Scope of the Declaration

The declaration applies to the use of the product in typical home, medical and commercial environments. The product should not be used in special electromagnetic environments.

Limitations of the declaration: The manufacturer is not responsible for damage or deterioration in product performance caused by improper use or use in special conditions inconsistent with the instructions for use.

TECHNICAL SPECIFICATION

Net weight: 18.0 kg

Maximum load: 150 kg

Seat depth: 57.9 cm

Seat width: 65 cm

Toilet opening width: 22.9 cm

Armrest height: 22.6 cm

Minimum height: 12 cm

Maximum lift – front: 21 cm

Maximum lift – rear: 37 cm

Maximum inclination angle: 15°

Minimum inclination angle: 0°

Water resistance protection rating: IPX4

Toilet opening length: 31.8 cm

Specification of the Wired Hand Controller

Rated input voltage: 100–240V, 50/60Hz, 2.0A

Rated output voltage: 29.0V, 1.5A

Specification of the Battery-Powered Hand Controller

Rated input voltage: DC 16.8V 2.2Ah

Rated output voltage: 29.0V 1.5A

Built-in battery: 14.8V 2.2Ah

Material composition:

powder-coated steel, stainless steel, PP/POM/nylon plastics, rubber, EVA, zinc alloy components.

WARNING! INDIVIDUAL PRODUCT DIMENSIONS MAY VARY WITHIN +/-2CM

This symbol indicates the maximum user weight.

PRODUCT COMPONENTS

- 1.Wired hand controller
- 2.Wired mains power adapter (A)
- 3.Toilet seat
- 4.Toilet lid and seat mounting holes
- 5.Toilet base
- 6.Battery-powered hand controller
- 7.Toilet base mounting holes
- 8.Battery-powered controller charger (B)
- 9.Padded armrests
- 10.Scissor mechanism

PARTS LIST

- A. Main body ×1
- B. Padded armrests ×2
- C. Wired hand controller ×1
- D. Wired mains power adapter (A) ×1
- E. Battery-powered hand controller ×1
- F. Battery-powered controller charger (B) ×1
- G. White hooks ×2

Note

- The wired mains power adapter (A) must be used only to connect and power the product.
- The wired hand controller requires direct connection to the control port of the raised toilet seat.
- The battery-powered controller charger (B) is used only for charging the battery-powered hand controller.
- These devices must not be used to power other products.

LIST OF INSTALLATION COMPONENTS: Expansion screw sets ×2, M8 self-locking nuts ×2, M8 open-end wrench ×1, Phillips screwdriver ×1, white nut covers ×2, white sleeves ×4

Note: Before installation, make sure that all components are complete and undamaged.

REQUIREMENTS FOR INSTALLATION ON THE TOILET

a. Load Capacity

The weight of the raised toilet seat is 18.0 kg. The maximum load of the raised toilet seat is 150 kg. Make sure that the total weight of the user and the raised toilet seat is within the permissible load capacity of the toilet.

b. Toilet Shape

The raised toilet seat is compatible only with round or elongated (oval) toilets.

c. Toilet Height

The raised toilet seat has no specific requirements regarding toilet height, which makes it suitable for toilets of various heights. The initial height of the raised toilet seat is 12 cm.

d. Toilet Dimensions

Before installation, measure the distance between the toilet mounting holes and the water tank, as well as the installation length (the distance from the mounting holes to the front edge of the toilet).

The toilet is compatible with the raised toilet seat if it meets one of the following conditions:

- Distance between the mounting holes and the tank: 2.5–4.3 cm

Installation length of the toilet: 42.4–55.9 cm


- Distance between the mounting holes and the tank: 4.3–6.4 cm

Installation length of the toilet: 40.6–53.8 cm

- Distance between the mounting holes and the tank: more than 6.4 cm

Installation length of the toilet: 38.6–50.8 cm

Toilets with a distance of less than 2.5 cm between the mounting holes and the tank are not compatible with the raised toilet seat.

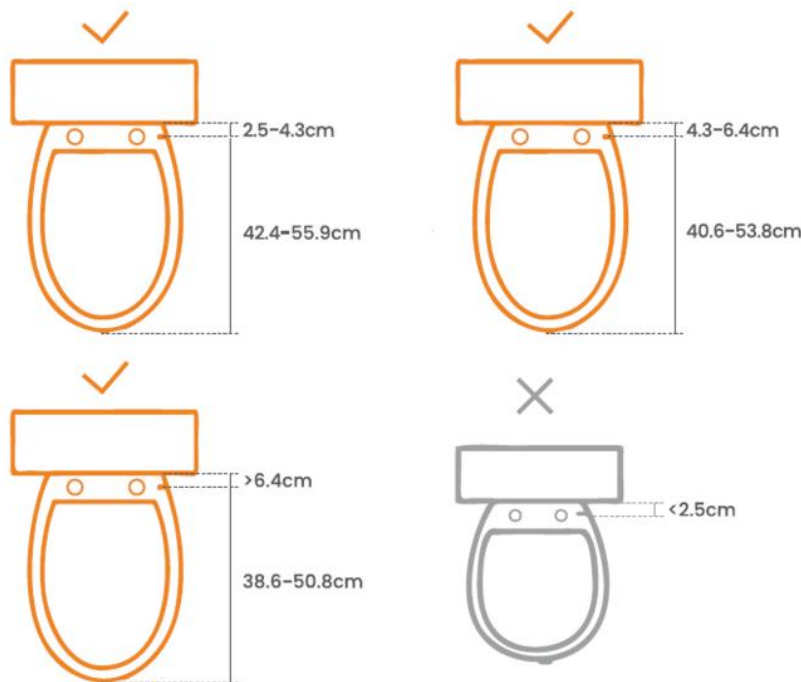REQUIREMENTS FOR INSTALLATION SPACE

a. Required Bathroom Entrance Width

The width of the raised toilet seat is 65 cm.

b. Required Installation Width

Make sure that the distance from the centre of the toilet to the walls on both sides is at least 37 cm.

c. Required Installation Length

The total length of the raised toilet seat may change slightly during the lifting process.

- In the lowest position, the length is 57.9 cm.

- In the highest position, the length is 62.4 cm.

d. Required Working Space

The distance between the raised toilet seat and elements such as the toilet paper holder should be at least 5 cm.


If it is not possible to maintain a clearance of 5 cm, move the toilet paper holder or other objects outside the range of movement of the raised toilet seat (e.g. outside the movement area of the seat or armrests) before installing the device.

e. Requirements Concerning the Power Cord and Socket

The length of the power cord of the raised toilet seat is 2 m.

Make sure that there is a safe and dry electrical socket nearby, within reach of the power cord.

It is recommended to use a waterproof socket with a protective cover.

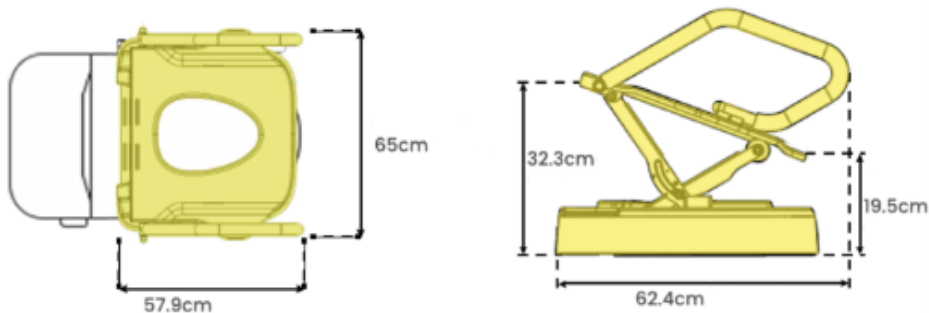INSTALLATION STEPS

Removing the Toilet Seat and Lid

Note: If the toilet seat and lid have already been removed and the toilet surface is clean and undamaged, proceed to the section: "Installing the raised toilet seat".

Prepare a screwdriver, wrench, protective gloves and a cloth.

Put on protective gloves to protect yourself from potential bacteria present on the toilet.

Prepare the working area: Make sure that there are no obstacles around the toilet seat and lid, which will make installation easier.

Locate the fastening elements: Check the rear part of the toilet seat and lid. Behind the seat there are two plastic or metal covers.

Open the covers: Use a screwdriver to pry open the covers to expose the screws underneath. The screws are secured with nuts.

Unscrew the screws: Hold the nuts by hand or with a wrench while using a screwdriver to unscrew the screws.

Remove the screws and nuts: After loosening the screws and nuts, put them aside in a visible place for later use.

Remove the toilet seat and lid: Carefully remove the toilet seat and the toilet lid.

Clean the toilet surface: Make sure that there is no dirt or remaining installation components left on the toilet surface. Check the toilet surface for damage. If damage is detected, repair it before further installation.


Installing the Raised Toilet Seat on the Toilet

1. Prepare the materials: Prepare the expansion screw sets and the Phillips screwdriver.


Each expansion screw set consists of three elements: plastic sleeve, metal washer, screw.

2. Install the plastic sleeves: Insert the plastic sleeves from the expansion screw set into the toilet mounting holes from above.

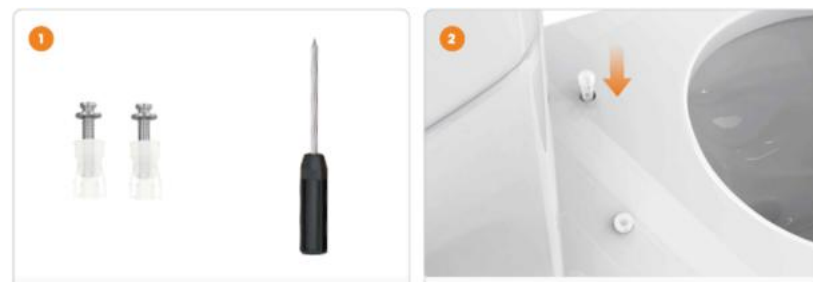3. Position the raised toilet seat: Place the raised toilet seat on the toilet, making sure that it is positioned stably and centrally.

4. Align the mounting holes: Make sure that the mounting holes in the base of the raised toilet seat are aligned with the toilet mounting holes. Select the appropriate installation position depending on the size of the toilet. The base of the raised toilet seat has three adjustable installation positions at the rear of the device.

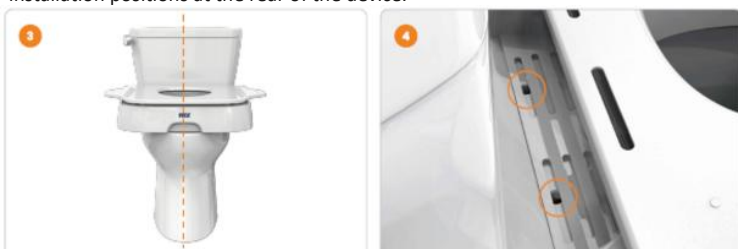5.Raise the raised toilet seat: After connecting the power supply, press and hold the “Up” button on the hand controller to raise the raised toilet seat to the highest position.

6.Place the metal washers: Place the metal washers from the expansion screw set on the mounting holes of the raised toilet seat base.

7.Fasten the screws: Pass the screws through the metal washers to attach the base of the raised toilet seat to the toilet mounting holes. Use the Phillips screwdriver to tighten the screws, making sure that the connection is stable and secure.

INSTALLING THE ARMRESTS

1. Prepare the Components

Prepare: 2 padded armrests, 2 M8 self-locking nuts, 4 white sleeves, 2 white nut covers, 1 M8 open-end wrench.

2. Install the Armrests

Place the armrest holes on the pivot axis located in the base of the device.

Secure the armrest by placing one M8 self-locking nut on the pivot axis.

Use the M8 open-end wrench to tighten the nut, then attach the white nut cover.

Perform the same steps to install the armrest on the other side of the device.

Note: If the user does not need to sit on the toilet from the side, the front part of the armrests may be directly attached to the seat. Using a Phillips screwdriver, loosen the screws located at the bottom of the front part of the armrests. Pass the screws through the seat and connect the armrest with the armrest on the other side, then tighten the screws with the screwdriver.

Before performing the above action, make sure that the seat is in the highest position.

Note: For hygienic reasons and for easier cleaning, it is recommended to reinstall the toilet seat and toilet lid. Use the original toilet screws and nuts to connect the seat and lid to the lifting mechanism of the raised toilet seat, then tighten the screws. Attach the white hook to the wall to hang the hand controller. It is recommended to install the hook in a convenient and easily accessible place for the user while using the toilet, which will increase convenience and safety of use.

HAND CONTROLLER CONNECTION METHODS

1. Connection Method for the Wired Mains Power Adapter (A)

Connect the wired mains power adapter (A) to the hand controller and the raised toilet seat according to the connector markings.

2. Connection Method for the Battery-Powered Controller Charger (B)

- Connect the battery-powered hand controller to the product.
- Connect the battery-powered hand controller to the battery-powered controller charger (B).

OPERATION TEST

Make sure that the product is stably installed on the toilet and that all screws have been properly tightened.

Connect the power cord to an electrical socket, making sure that the place of use is dry and safe. Check that the power cord is not at risk of being pressed by other objects.

Operation of the Hand Controller

The hand controller has two buttons:

- press and hold the “Up” button to raise the raised toilet seat,
- press and hold the “Down” button to lower the raised toilet seat.

Momentary Control Mode

The controller operates in momentary control mode, which means that the button works only while it is being pressed. When the button is released, the device stops. In the event of an unexpected situation during operation, release the button immediately to stop the device.

NOTE: INITIAL RESET

Before first use, reset the device. Press and hold the “Up” and “Down” buttons simultaneously and continue holding them. The green indicator light will start flashing until a sound signal is heard from the motor. This means that the system has entered reset mode.

After the reset has been completed correctly, the green indicator light will remain on, which means that the device is ready for operation. If the reset fails, the green indicator light will continue flashing or will go out.

Operation Test

After completing the reset, perform a test of raising and lowering the device. Make sure that the raised toilet seat (including the seat and armrests) has sufficient working space without obstacles such as walls, toilet paper holders, water tank or other equipment. Make sure that the user can conveniently operate and adjust the raised toilet seat.

WIRED HAND REMOTE CONTROL

“UP” button: Press and hold the button to raise the device continuously.

The device will stop automatically after reaching the highest position.

“DOWN” button: Press and hold the button to lower the device continuously.

The device will stop automatically after reaching the lowest position.

Indicator light: A green light indicates correct operation of the device. If the indicator light is red, the reset procedure must be performed.

BATTERY-POWERED HAND REMOTE CONTROL

Indicator light: A green light indicates correct operation of the device. If the indicator light is red, the reset procedure must be performed.

STOP button: Press the “STOP” button to quickly mute the sound signal.

“UP” button: Press and hold the button to raise the device continuously.

The device will stop automatically after reaching the highest position.

“DOWN” button: Press and hold the button to lower the device continuously. The device will stop automatically after reaching the lowest position.

Motor Zero Position Calibration / Reset Function

Press and hold the “UP” and “DOWN” buttons simultaneously to start the reset process.

Continue holding the buttons until the motor emits a sound signal.

Battery Warning

When the battery level of the remote control drops below 30%, the device will emit three short sound signals reminding the user that charging is required.

When the battery level drops below 10%, a continuous sound signal will be emitted approximately every 3 seconds, indicating a critically low battery level.

Charge the device as soon as possible to prevent excessive discharge during a longer period of non-use, which could make it impossible to recharge the battery.

TROUBLESHOOTING

Problem: After pressing the button, the seat does not raise or lower.

Check:

- Check whether the power socket is not damp and whether the power cord is properly connected.
- Check whether the indicator light of the hand controller is on.

Solution: If the device does not start or there are problems with operation, make sure that the hand controller is powered before performing the reset.

Reset procedure:

Press and hold the “Up” and “Down” buttons simultaneously until a continuous sound signal is heard.

The sound signal means that the system has entered reset mode.

After entering reset mode, the green indicator light will start flashing.

After the reset has been completed correctly, the green indicator light will remain steadily on during normal operation.

If an error occurs, the red indicator light will start flashing, indicating a problem.

If the reset fails, contact the seller.

USE

1. **Raising the Seat:** Use the hand remote control to raise the raised toilet seat to the appropriate height to make sitting down easier.

2. **Sitting Down Slowly:** The user should sit down slowly and calmly in the centre of the seat. Avoid sitting down suddenly in order not to cause malfunction of the device.

3. **Holding the Armrest and Remote Control:** The user should hold the armrest steadily with one hand and hold the hand remote control with the other hand, resting it on the leg or armrest. Then press and hold the “DOWN” button on the hand remote control.

4. **Lowering to a Comfortable Position:** After lowering the device to a comfortable position, release the button on the remote control and use the toilet.

5. **Raising After Use:** After using the toilet, press and hold the “UP” button on the hand remote control until the device rises to a comfortable height allowing the user to stand up.

6. **Standing Up Safely:** The user should grasp the armrests with both hands and use their support to stand up safely.

NOTE: In the event of a “serious incident” related to the device, which directly or indirectly led, could have led or may lead to any of the following events:

- a) death of the patient, user or another person, or
- b) temporary or permanent deterioration of the health condition of the patient, user or another person, or
- c) a serious public health threat, the above-mentioned “serious incident” must be reported to the manufacturer and to the competent authority of the Member State in which the user or patient resides. In the case of the United Kingdom, the competent authority is the Medicines and Healthcare products Regulatory Agency (MHRA). **NOTE:** In the event of pain, allergic reactions or other disturbing symptoms related to the use of the medical device that are unclear to the user, consult a healthcare professional. **NOTE:** It is forbidden to use the product in a manner other than in accordance with its intended use.

MAINTENANCE AND CLEANING

- It is recommended to check the device regularly to ensure that it works properly. Before each use, check the condition of the product, cleanliness, condition of screws, caps, etc.
- Before cleaning the product, make sure that the power supply has been disconnected.
- Use a soft cloth and a non-corrosive cleaning agent to clean the product. The surface should be wiped gently to avoid scratches or damage to the product.
- Do not use cleaning agents containing alcohol to clean the internal moving parts and the motor. Direct washing may cause loss of lubricant and shorten the service life of the product. Avoid moisture entering the internal moving parts and the motor.
- After each use by another person, the product must be thoroughly cleaned and disinfected.
- Avoid contact with sharp objects to prevent paint chipping. **Note:** If visible damage, unusual sounds or malfunction are noticed, stop using the product immediately.

STORAGE: Do not expose to direct sunlight. Storage temperature and humidity: -10°C to +60°C; relative humidity from 0% RH to 75% RH.

DISPOSAL: The used product must be disposed of in accordance with local regulations concerning waste electrical and electronic equipment (WEEE). The product contains electrical components and a battery, therefore it must not be disposed of together with municipal waste. The battery and electrical parts should be taken to an appropriate collection point for electrical waste or recycling. Batteries must not be thrown into fire or exposed to mechanical damage, as this may cause fire, leakage of hazardous substances or explosion.

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Model:
Serial number:

Point of sale stamp

and a legible signature of the Seller:

1. ANTAR MEDICAL Sp. z o.o. ,03-068 Warsaw, ul. Zawisłańska 43, hereby grants a 12-month warranty on the product, starting from the date of purchase by the Buyer.
2. During the warranty period ANTAR undertakes to carry out the repairs necessary to restore proper functioning of the product, within 14 days from the date of receipt of the product by the service centre, together with the proof of purchase and the warranty card.
3. ANTAR reserves the right to replace the product if the repair is deemed uneconomic. Repair or replacement of faulty products does not extend the warranty period.
4. If the product is found to be malfunctioning, please contact your retailer/distributor or (if it is impossible) send it to ANTAR (as specified in point 1).
5. Before use, read the instructions carefully and follow them. If the product is used improperly and contrary to recommendations, the warranty will be irrevocably lost.
6. The warranty does not cover damage caused by external mechanical force, load exceeding the maximum allowed weight of the user, as described in the manual, as well as contact of the product with alcohol, greasy substances or petrol. The warranty does not cover parts which are subject to normal wear or tear during use, such as the bearing, rubber parts and others.
7. The product delivered to the service in connection with the complaint cannot cause an epidemic hazard.
8. THIS WARRANTY IS ONLY VALID TOGETHER WITH THE PROOF OF PURCHASE (BILL, RECEIPT, INVOICE). IN CASE OF A COMPLAINT, A PROOF OF PURCHASE MUST BE ATTACHED.
9. Complaints submitted without the proof of purchase and a warranty card with the product serial number will not be considered.
10. The warranty does not cover parts and components subject to natural wear and tear, such as wheels, tires, rubber pads, fuses, etc., the wear of which results from normal use of the product. The warranty period for batteries is 6 months.

ANTAR MEDICAL Sp. z o.o.
03-068 Warsaw, ul. Zawisłańska 43
Tel. 22 518 36 00, Fax 22 518 36 30
www.antar.net

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